ZOLID-SHT

ZOLID-SHT preshade

ZOLID-SHT multilayer

DE	Gebrauchsanweisung	3 - 21
EN	Instruction Manual	22 - 40
FR	Instructions d'utilisation	41 - 60
IT	Istruzioni per l'uso	61 - 79
ES	Instrucciones de uso	80 - 100
PT	Instruções de uso	101 - 119

- Original Gebrauchsanweisung -

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- Translation of the original Operating Instructions -

Table of Contents

Explanation of Symbols	23
General Safety Instructions	25
Suitable Personnel	25
Properties	25
Warranty/Exclusion of Liability	27
Testing the Blanks	28
Mounting to the Worktable	28
Application	28
Fixation	38
Post-processing at the dentist's	39
Follow-up examination at the dentist's	39
Environmental Protection	40
Accessories	40
Download information	40

Explanation of Symbols

Warning indications

Warning indications in the text are marked with a triangle and boxed.

In case of hazards through electricity, the exclamation mark in the warning triangle is substituted by a lightning bolt.

Signal words at the beginning of a warning indication specify the type and severity of the consequences, if the measures to avert the hazard are not adhered to.

- _ **NOTE** means that property damage can occur.
- _ **CAUTION** means that light to fairly serious personal injury can occur.
- _ **WARNING** means that serious personal injury can occur.
- _ **DANGER** means that life-threatening personal injury can occur.

Important information

Important information that do not lead to hazards for humans or to property damage are marked with the icon aside and are boxed.

Other symbols in the Manual

Symbol	Meaning
▷	Item of an operation description
—	Item of a list
•	Subitem of an operation description or a list
[3]	Numbers in square brackets refer to position numbers in graphics/figures

Other symbols on the product

Symbol	Meaning
	Catalog number
	Batch code
	Manufacturer
	Consult instructions for use
	Used by YYYY-MM-DD or YYYY-MM
Rx only	Symbol that may be used in place of statement “CAUTION”. US Federal law restricts this device to sale by or on the order of a dentist.
	Serial number

General Safety Instructions

CAUTION:

Possible health impairment from zirconium-oxide dust!

- ▷ When processing, wear personal protective equipment (dust protection mask, safety glasses/goggles, ...).
- ▷ Do not breathe dust/fume/gas/mist/vapours/spray.
- ▷ Avoid release to the environment.

Suitable Personnel

NOTE:

The product may only be used by trained dental technicians.

Properties

Intended use

Zirconium-oxide blanks for permanent and removable dental prosthetics.

Product Description

The blanks are made of zirconium oxide (ZrO_2) for type II, class 4 dental applications in accordance with DIN EN ISO 6872.

They are used for manufacturing permanent and removable prosthetic restorations (e.g., crowns and three-unit bridges to the molar region) using CNC milling machines (e.g., Straumann milling machines).

After completion of the specified end-sintering, the blanks meet the requirements of DIN EN ISO 6872.

Safety data sheet/Declaration of conformity

The data sheet on safety and the declaration of conformity are available on request.

Technical data

	Unit	ZOLID-SHT ZOLID-SHT preshade ZOLID-SHT multilayer
Bending strength (3-point)	MPa	700 ± 150
Bending strength (4-point)	MPa	600 ± 150
E-module	GPa	≥ 200
Thermal expansion coefficient (CTE) (25 - 500°C)	10 ⁻⁶ /K	10.1 ± 0.5
Chemical solubility	µg/cm ²	< 100
Vickers hardness	HV10	1300 ± 200

The delivery scope can vary country-specific.

Chemical composition

Oxide	Mass percentage		
		ZOLID-SHT	ZOLID-SHT
	ZOLID-SHT	preshade	multilayer
$ZrO_2 + HfO_2 + Y_2O_3$	≥ 99.0	≥ 99.0	≥ 99.0
Y_2O_3	9.15 - 9.55	9.15 - 9.55	8.5 - 9.5
HfO_2	≤ 5	≤ 5	≤ 5
Al_2O_3	≤ 0.5	≤ 0.5	≤ 0.5
Other oxides	≤ 1	≤ 1	≤ 1

Shelf life of the blanks

When appropriately stored, the blanks can be used 5 years after the manufacturing date.

Storage

Store the blanks in the original packaging in a dry location.

Warranty/Exclusion of Liability

Application-technical recommendations, whether given orally, in writing or in the course of practical training, are guidelines. Our products are subject to continuous further development. Subject to changes in handling and composition.

Testing the Blanks

Upon receipt, it is imperative to visually check the proper condition of the blanks. After using a (transport-) damaged blank, complaint claims are no longer possible.

Mounting to the Worktable

The mounting of the blanks is described in the Operating Instructions of the respective milling machine.

Application

Indication ranges

- _ Anatomically reduced crown framework and three-unit bridge frames in the anterior and posterior tooth range.
- _ Monolithic (fully anatomical) crowns and up to three-unit bridges in the anterior and posterior tooth range.
- _ Anatomically reduced bridge frames with a maximum of one intermediate link in the anterior and posterior tooth range.
- _ Monolithic bridges with a maximum of one intermediate link in the anterior and posterior tooth range.

Contraindications

- _ Insufficient tooth-structure availability
- _ Insufficient preparation results
- _ Insufficient oral hygiene
- _ More than one bridge unit
- _ Known incompatibilities with respect to the components
- _ Heavily discoloured hard tooth structure
- _ provisional insertion

Material-specific frame parameters

The following material-specific frame parameters must be complied with when fabricating zirconium oxide frames in high-density sintered condition:

Minimum frame thickness in mm	Connector cross-section in mm ²		Maximum number of connected bridge units	
	Anterior	Posterior	Anterior	Posterior
0.5	≥ 12	≥ 12	1	1

Detailed minimum wall thicknesses and connector cross-sections

Indication	Pattern	Total number of units	Number of connected bridge units	Wall thickness in mm		Connector cross-section in mm ²
				incisal/occlusal	Circular	
Primary components/Double crowns	–	1	–	0.7	0.5	–
Single crown	–	1	–	0.5	0.5	–
Front-tooth bridge	OXO	3	1	0.5	0.5	≥ 12
Posterior bridge	OXO	3	1	0.7	0.5	≥ 12


○ Abutment crown

✕ Bridge unit

Crown border parameters in the CAD software

Pay attention to sufficient crown border thicknesses Recommended values:

- 1. Crown border thickness (Horizontal): 0.2 mm
- 2. Angled: 0.3 mm
- 3. Angle: 60°
- 4. Height (Vertical): 0 mm

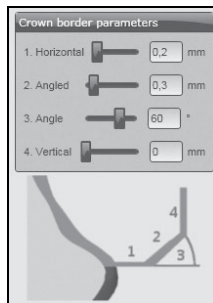For tangential preparations, increase the “crown border thickness” and “Angled” parameters, e.g.:



- _ 1. Crown border thickness (horizontal): 0.25 - 0.3 mm
- _ 2. Angled: 0.35 - 0.4 mm

Observe zirconium-oxide-efficient design forms and do not create pointed edges.

Positioning in the CAM software

- ▷ Set adequate holding ridges with sufficient thickness (at least 3 holding ridges for single crowns; thickness > 1.7 mm).
- ▷ Position holding ridges approx. 1 - 2 mm away from the crown margin, especially for long crown margins. This provides support to the crown margin.
- ▷ Position the crown in the blank in such a manner that the clearance from the crown margin to the blank surface is at least 1 mm. This ensures that the crown has no contact to the sinter skin of the blank.

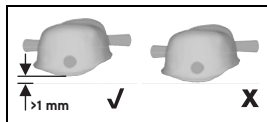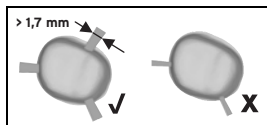Procedure when fabricating dental framework

For maximum fitting accuracy, the specific expansion/shrinkage value is given on the blanks.

When machining the blanks in milling machines, it is recommended to only use cutters with sufficient cutting edge quality.

- ▷ Depending on requirement in the respective CAM software, enter the appropriate value.

In case of uncertainties, contact your CAM supplier.

- ▷ Mill the blanks with the milling machine applying the milling strategies provided for the material.

Removal from the blank

Utmost care must be taken when removing the frames from the blank:

- ▷ Separate the frames by grinding them out using a turbine grinder (without water-cooling) and a pointed diamond bit (do not use cutting discs).

For trimming and dressing the connectors and smoothing out the frames, we recommend using a polishing set.

- ▷ Clean the frames from residual zirconia dust using a brush or compressed air.

Colouration of the white zirconium oxide

- ▷ Depending on requirement, partially or completely stain the white zirconium oxide frames prior to the sintering with dye solution.

Colouration of the pre-stained monochrome zirconium oxide

These blanks are not stained with the dye solutions, but monolithically stained upon successful final sintering with stain or glazing material or used as a frame for veneering.

Blank Shade	Tooth shade							
	Bleach	A1	A2	A3	A3,5	B1	B2	B3
ZOLID-SHT pre-shade Bleach	●	-	-	-	-	-	-	-
ZOLID-SHT pre-shade A light	-	●	●	○	-	-	-	-
ZOLID-SHT pre-shade A medium	-	-	●	●	○	-	-	-
ZOLID-SHT pre-shade B light	-	-	-	-	-	●	●	○

- ideal for this shade
- shade can be achieved
- not suitable for this shade

Further shades are used in the same manner

Colouration of the pre-stained polychrome zirconium oxide

These blanks are not stained with the dye solutions, but monolithically stained upon successful final sintering with stain or glazing material or used as a frame for veneering.

Shade selection and placement of the restoration in the blank

The tooth shade of the blank covers off the imprinted area (e. g. 0/A1). When a higher blank is selected, the respective lighter or darker shade can be achieved more significantly by vertically moving in the nesting. The closer the work is positioned to the imprinted side, the larger the lighter shade portion will be.

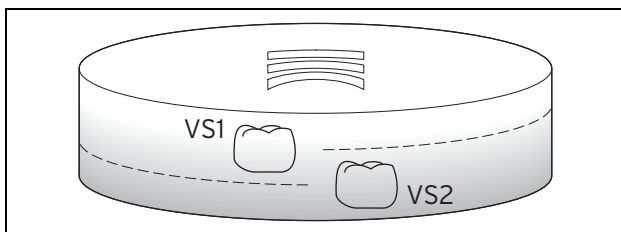
Fig. 2 Placement of the restoration in the blank

VS1 Tooth shade 1 (lighter)

VS2 Tooth shade 2 (darker)

By applying the stain, further tooth shades can be achieved.

Final sintering

A zirconium sintering furnace or another high-quality and suitable furnace is recommended for the sintering (please also refer to the operating instructions of the respective sintering furnace).

Highly translucent zirconium oxide and opaque zirconium oxide can be sintered together in one sintering process.

**NOTE:**

When highly translucent zirconium oxide and opaque zirconium oxide are sintered together in the same sintering bowl and with the same sintering pearls, then this can lead to changes in the translucency. Therefore, it is advisable for highly translucent zirconium oxide and opaque zirconium oxide, to use individual sintering pearls in separate sintering bowls.

- ▷ Place the frames for sintering into the sintering bowl filled with sintering pearls.
- ▷ To ensure proper support, place the frames with slight pressure onto the sintering pearls.

**NOTE:**

▷ Pay attention that no pearls are jammed in the interdental spaces or located in the crown cavities!

- ▷ Perform the final sintering of the frames according to the following program:
 - Heat-up phase: Room temperature to final temperature 1450 °C; heating rate 5 - 10 K/min
 - Dwell time at final temperature: 2 hrs
 - Cooling phase: Final temperature 1450 °C to room temperature (at least < 200 °C); approx. 5 K/min (approx. 5 hours)

Frames for veneering

Post-processing

If required, the zirconium oxide frames can be reworked after the final sintering without pressure, using a water-cooled lab turbine and suitable diamond grinding points (recommended grain approx. 40 μm).

Frame preparation for veneering

After checking the fitting accuracy, contact points and occlusion:

- ▷ Clean the interior crown surfaces by means of corundum blasting (aluminium oxide, 110 μm , pressure ≤ 2 bar) and subsequent steaming off.

Veneer ceramics

For information on the thermal-expansion coefficient of the blanks, please refer to the Technical data (see page 26).

Monolithic dentures

Post-processing

After the final sintering has taken place:

- ▷ If required, the zirconium oxide frames can be reworked using a water-cooled lab turbine and suitable diamond grinding points (recommended grain approx. 40 μm).
- ▷ Polish the contact surfaces to the antagonist and to the neighbouring teeth to a high gloss when in high-density sintered condition.

Preparation for staining and glazing

After checking the fitting accuracy, contact points and occlusion:

- ▷ Clean the interior crown surfaces by means of corundum blasting (aluminium oxide, 110 μm , pressure ≤ 2 bar) and subsequent ultrasonic treatment and/or steaming off.

The exterior crown surfaces are not corundum blasted.

For bridges with massive intermediate units, we recommend step-by-step long-term cooling to 500 °C.

Staining and glazing

- ▷ Stain and glaze the units afterwards.

Fixation

Due to the high translucency of this zirconium oxide, the selection of cement shade and possible discolourations of the tooth stump are to be observed, as these may affect the final colour result.

Conventional cementation:

Due to the high strength and stability of the frames, conventional fixation with zinc-oxide phosphate cement or glass ionomer cement (e.g., GC Fuji Plus) is possible in most cases.

► When applying conventional cementation, observe sufficient retention and a minimum abutment height of 3 mm!

Adhesive fixation:

For adhesive fixation, we recommend the bonding composite Panavia® 21 or Panavia® F 2.0, owing to the excellent bond strength with the zirconium-oxide frame material.

Provisional cementation is not recommended, as the frames can be subject to damage upon removal.

Post-processing at the dentist's

When grinding is required by the dentist while installing the work into the patient's mouth, the units must be well re-polished to a high gloss. For grinding, we recommend using only diamond grinding points (recommended grain approx. 40 μm). For polishing, we recommend using diamond polishing bits.

NOTE:

When the work is not sufficiently polished, the antagonist may become damaged due to abrasion!

Follow-up examination at the dentist's

It is urgently recommended to check the monolithic dentures in the patient's mouth once a year. In this, the remaining dentition, the antagonists and the soft tissue are also to be inspected. Corrective measures are to be undertaken as required. Here, it is also important that the objects are polished afterwards to a high gloss.

Environmental Protection

Packaging

In terms of packaging, the manufacturer participates in country-specific recycling systems, which ensure optimal recycling.

All packaging materials used are environmentally-friendly and recyclable.

Disposal of the blanks

Larger quantities may not be disposed of together with household waste. Smaller quantities may be disposed of together with household waste. Do not dispose of into sewer systems. Dispose of non-cleanable packaging properly.

Always dispose of according to official regulations.

Accessories

For more information on product-specific accessories, please refer to www.amanngirrbach.com or www.straumann.com.

Download information

Further instructions/manuals can be downloaded under ifu.straumann.com.

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Amann Girrbach AG
Institut Straumann AG