# ST.330.038

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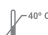
English  Pilot Drill

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Simbologia Symbology Simbología Simbologia Symbolik Symboles	Descrição / Description / Descripción / Descrizione / Beschreibung / Description
	Tamanho do produto / Product size / Tamaño del producto / Dimensioni del prodotto / Produktmaße / Taille du produit
	Código do produto / Product code / Código del producto / Codice del prodotto / Produktcode / Code du produit
	Número do lote / Batch Number / Número de Partida / Numero di lotto / Chargennummer / Numéro de lot
	Data de fabricação / Date of manufacture / Fecha de fabricación / Data di fabbricazione / Herstellungsdatum / Date de fabrication
	Prazo de validade / Shelf life / Plazo de validad / Durata di conservazione / Mindesthaltbarkeit / Date de péremption
	Quantidade / Quantity / Cantidad / Quantità / Qualität / Quantité
	Material de fabricação do produto / Material used / Material de fabricación del producto / Materiale di fabbricazione del prodotto / Herstellungsmaterial des Produkts / Matériau de base du produit
	Não estéril / Non-sterile / No estéril / Non sterile / Nicht steril / Non stérile

Simbologia Symbology Simbología Simbologia Symbolik Symboles	Descrição / Description / Descripción / Descrizione / Beschreibung / Description
	Consulte as instruções de utilização / Refer to instructions for use / Consulte las instrucciones de utilización / Consultare le istruzioni per l'uso / Siehe Bedienungsanleitung / Consulter les instructions d'utilisation
	Limite superior de temperatura / Upper limit of temperature / Limite superior de temperatura / Limite di temperatura massima / Maximaltemperatur / Limite de température supérieure
	Conservar seco / Keep dry / Consérvelo seco / Conservare in luogo asciutto / Trodden lagern / Garder au sec
	Manter afastado da luz solar / Keep protected from sunlight / Manténgalo lejos de la luz solar / Tenere al riparo dalla luce solare / Vor Sonneneinstrahlung schützen / Conserver à l'abri de la lumière solaire
	Não utilizar se a embalagem estiver danificada / Do not use if package is damaged / No lo utilice si el envoltorio está dañado / Non utilizzare se la confezione originale è danneggiata / Bei beschädigter Verpackung Produkt nicht verwenden / Ne pas utiliser si l'emballage a été abîmé

Simbologia Symbology Simbología Simbologia Symbolik Symboles	Descrição / Description / Descripción / Descrizione / Beschreibung / Description
	Marcação CE para comercialização na Comunidade Europeia / CE Mark for European Community market / Marca CE para comercialización en la Comunidad Europea / Marchio CE per la commercializzazione in Europa / CE Kennzeichen zum Vertrieb in der Europäischen Gemeinschaft / Marquage CE pour une commercialisation en Europe
	Notificação exigida pelo FDA para comercialização nos Estados Unidos / Notification required by FDA for United States market / Notificación exigida por el FDA para comercialización en los Estados Unidos / Notifica necessaria per la commercializzazione negli Stati Uniti / Vom FDA geforderte Kennzeichnung für den Vertrieb in den USA / Notification exigée par FDA pour une commercialisation aux États-Unis
	Representante na Comunidade Europeia / Representative in the European Community / Representante en la Comunidad Europea / Rappresentante nella Comunità europea / Vertreter in der Europäischen Gemeinschaft / Représentant dans la Communauté Européenne
	Fabricante / Manufacturer / Fabricante / Produttore / Hersteller / Fabricant

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This device is intended for a specialized procedure, which should be performed by professionals qualified in Dental Implants. For optimum results, use the product knowing the appropriate techniques. Always apply them under appropriate conditions, also in an operating room atmosphere.

## DESCRIPTION

The Pilot Drill is manufactured in surgical stainless steel. In its active end there is a reduced guide diameter, for positioning during drilling and, cutting edges in its nominal diameter. In the other end there is fitting for Contra-Angle in accordance to standard ISO 1797-1 - Dental rotary instruments - Shank. It is available in two presentations: 2/3 Pilot Drill and CM Pilot Drill. The laser markings determine: for the 2/3 Pilot Drill - the final depth of the drilling and, for the CM Pilot Drill - 3 installation levels of the CM Implant, at level, 1 or 2mm below the bone crest. The Pilot Drill receives the application of a dark color thin carbon film on the cutting surface. This film has as main purposes: • reduction of the heat generated during osteotomy • friction reduction between the Drill and the bone • resistance increase to wear • resistance increase to oxidation.

## APPLICATIONS

This product is a surgical instrument used for drilling the bone tissue, with the purpose of:  
 - Intermediary drilling, in the case of the 2/3 Pilot Drill, with the function of increasing the 2.0mm drilling diameter and serving as guide for the following 3.0mm Drill entrance or,  
 - Final drilling, in the case of the CM Pilot Drills, with the function of preparing the bone crest for the settlement of the cervical third of the Morse Taper Implant chosen.

## CONTRAINDICATIONS

This product does not present contraindications, once they are properly used for the indicated purposes.

## HANDLING

Fit the Drill in the Contra-Angle. Place the Drill in the specific rotation. Perform the drilling. During drilling, the pressure should not be excessive and, an abundant irrigation, either manual or combined, using the handpiece pump is necessary. Note 1: In order to minimize the surgical trauma, it is important to respect the Drill sequence for the success of the Implant osseointegration. Refer to Neodent catalog or the website [www.neodent.com.br](http://www.neodent.com.br). Note 2: Do not interrupt the handpiece rotation with the Drill within the surgical cavity, because this can make difficult the removal or cause the Drill break.

## SANITATION

The surgical instruments should be correctly sanitized after each use. Do not leave the product in contact with humidity more than the time necessary for the sanitation. Proceed as follows: 1<sup>st</sup> step: Fully immerse the part in enzymatic detergent (diluted according to the manufacturer). 2<sup>nd</sup> step: Wash in ultrasonic cleaner for approximately 10 to 15 minutes. 3<sup>rd</sup> step: Rinse with distilled water in abundance, until completely remove the wastes from the solution. Nylon brushes use is recommended. 4<sup>th</sup> step: Dry with a clean and dry cloth or with compressed air. 5<sup>th</sup> step: Perform a visual inspection, observing if there are failures in the cleaning process. If there are still wastes, the part should be once again immersed in detergent - 1<sup>st</sup> step - and, if necessary, the cleaning should be performed with the aid of a nylon brush. Repeat the rinsing and drying sequence. 6<sup>th</sup> step: If applicable, place the product in an appropriate case duly sanitized and dry. 7<sup>th</sup> step: Select the packing in accordance to the sterilization process. Use preferably surgical grade paper self-sealing packing with laminated film.

## PRESENTATION AND STERILIZATION

This product is unitarily provided, conditioned in blister (surgical grade film and paper type) packing. It is reusable, delivered non-sterile and should, therefore, be sterilized before use. Sterilize the product on the day before or on the day of the procedure. We preferably recommend following the sterilization method by steam autoclaving and, yet, the parameters and procedures established in standard ISO 17665-1:2006 - Sterilization of health care products. Moist heat. Requirements for the development, validation and routine control of a sterilization process for medical devices. Suitable parameters: 121°C, 30 min, at 1 atm. Alternatively, can proceed according to the instruction manual of the manufacturer of the autoclave. Sterilization validity: 7 to 15 days, provided that the product is conditioned in clean, dry and protected from the sunlight environment.

## PRECAUTIONS

• Do not use the product if the packing is damaged. • This product should be sterile. • The non-replacement of the Drills, as recommended by the manufacturer, can generate improper bone heating, compromising the osseointegration success. • Do not regrind. • Do not use solutions with ammonia, acids or hydrogen peroxide, because they can damage the Drill. • Do not use anticrust products, avoid using enzymatic liquid in concentrations above 10% and dry parts that still have cleaning solution wastes, because these procedures cause oxidation. • The use of

steel brushes is inadvisable. We recommend the use of nylon brushes. • The improper planning can jeopardize the performance of the Implant/prosthesis set causing failures to the system, such as loss or break of the Implant, loosening or fracture of the prosthetic screws. • Better results are obtained with the use of the Neodent products sequence. • The use of instruments and/or prosthetic abutments of other manufacturers do not guarantee the perfect function of the Neodent Implants System and exempts any guarantee of the product.

## ADVERSE EFFECTS

Adverse effects will only occur if the choice and the use of instrument are improper.

## POST-OPERATIVE PRECAUTIONS AND MAINTENANCE

Instruct the patient as to the need of a professional medical monitoring after the surgery and to obey the guidelines regarding the precautions, hygiene and prescription of drugs. These guidelines are the responsibility of the professional in charge.

## STORAGE CONDITIONS

This product should be stored in a clean and dry location, in a maximum temperature of 40°C and protected from direct sunlight.

## DISPOSAL OF MATERIAL

Every product and consumable used during the surgery for the installation of dental implants may endanger the health of those who handle them after use. Before discarding them into the environment, it is recommended to take a look at the current legislation and adhere to it.

## DATE OF EXPIRATION

Written on the label.

## LIFE CYCLE

This product has a life cycle of up to 20 drillings for bone qualities I and II, and up to 30 drillings for bone qualities III and IV, since its use conditions are respected, including adequate irrigation and rotation. Apart of the life cycle period mentioned, the professional should check the sharpening conditions of the Drill. Drills with improper cutting conditions should be disposed. Note: Classification of the bone quality according to Lekholm and Zarb (1985).

Federal (USA) law restricts this device to sale by or on the order of a licensed dentist or physician.

Not all products are available in all countries. Please, contact the authorized distributor.

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